

# Plain and Simple ICU Perupa

- Perupa can be installed with the cage that is in use.
- Well designed for reduce the routine maintenances.
- Obtain ICU as low cost.
- Available for use as HBO(Hyperbaric Oxygen Therapy).




- 1 Control :** Controlled Temperature and Oxygen concentration in the care room by easy operation buttons.
- 2 Safety :** Ventilator opens automatically when shutdown the power.
- 3 Lighting :** LED light is equipped for the care room, It has brightness with long life.
- 4 Dehumidification :** AC condensation drain line is equipped.

Specification of this product would be changed without any notice.


## Dimension

This dimension is Perupa ICU with B1/2 ST cage by TOKYO MENIX. This cage is Option.


\*Dimension (especially depth) is varied on the size of the optional cage.

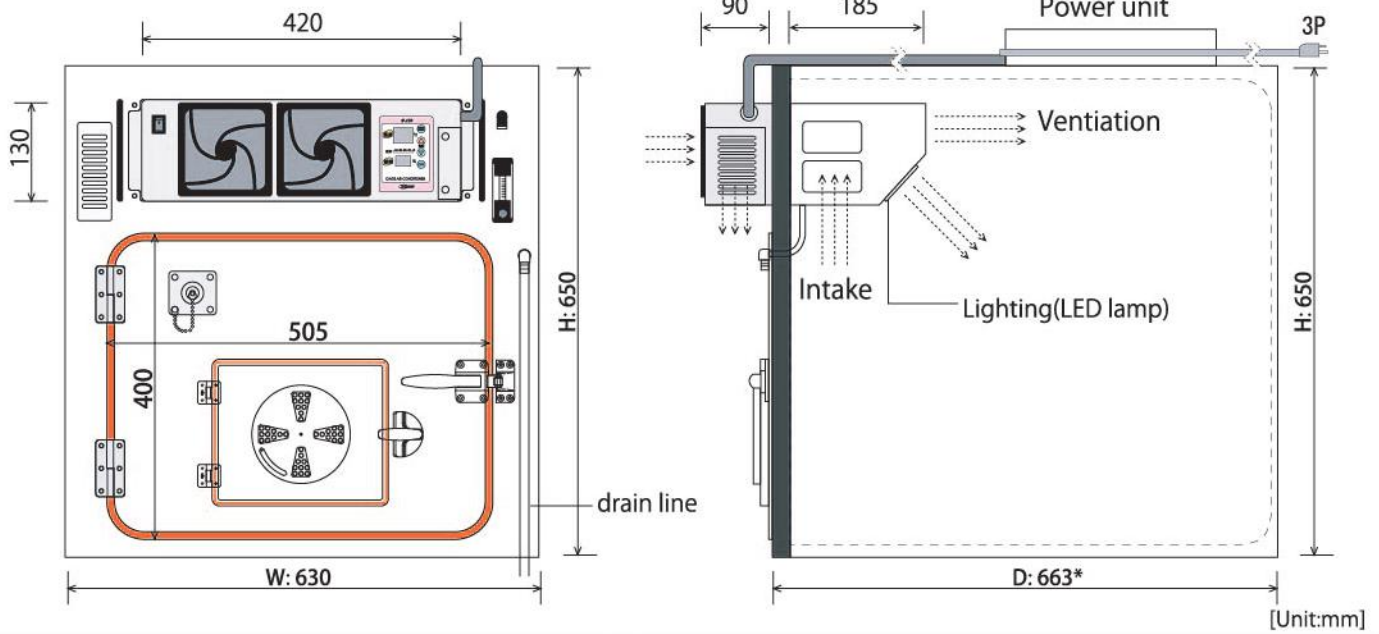

## Operation panel

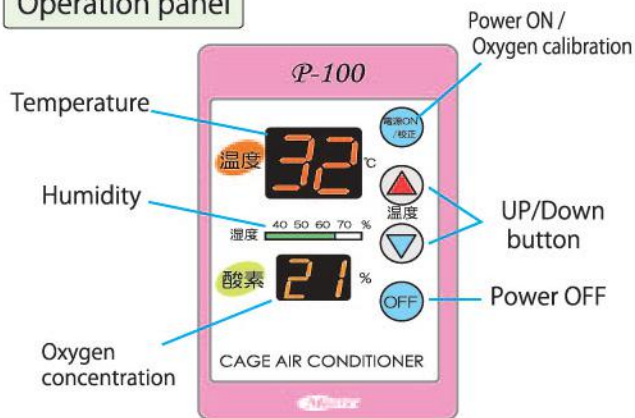

## Option

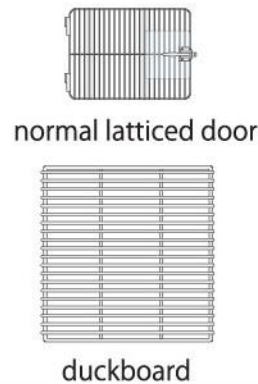

Control temperature range: between 17C to 35C	The range of temperature in the care room is varied on the condition of outer temperature or the animal that is in the care room.
Dehumidification range: between 40% to 60%	It's automatic function to control humidity between 40% to 60%RH when approx between 1kg to 5kgs animal is in the care room.
Control oxygen concentration range: between 21% to 60%	Electric oxygen controller is watching the condition constantly to adjust the oxygen concentration level.
Emergency ventilation system (Safety function)	Emergency vent opens automatically when shutdown the power supply. Perupa provides safety on sudden blackout.
Drain line	Perupa has AC condensation drain line to relieve water.
Quiet	Cooling fans are controlled by the condition of temperature in the care room.
Save energy	Reduce running cost, this is using thermoelectric cooling device for cooling the temperature (On cooling: 220W, Heating: 180W, Retaining temp: 80W)
Small oxygen flowmeter	Equipped small oxygen flowmeter to adjust the volume of oxygen into the care room.
Access door with ventilation	Equipped vent holes on access door. That would lead the fresh air into the care room

**Menix** TOKYO MENIX CO.,LTD.

2-1442-8, Kotesashi-cho, Tokorozawa-shi, Saitama-ken, 359-1141, Japan

TEL: +81-4-2923-8841 FAX: +81-4-2923-8864

http://www.t-menix.com e-mail: inter@t-menix.com